

UVT-B-AR, DNA/RNA UV-cleaner box

DESCRIPTION

DNA/RNA UV-cleaner box **UVT-B-AR** is designed for clean operations with DNA samples. UV-Cabinet provide protection against contamination.

Model is a bench-top type, made of metal framework, working surface made of stainless steel.

UV-cleaner boxes are equipped with an open UV lamp installed in the upper hood. UV-radiation from the open lamps disinfects the working area inactivating DNA/RNA fragments during 15–30 min of exposure. A digital timer controls duration of the direct UV irradiation. A daylight lamp provides proper illumination of the working surface.

UV-cleaner boxes are equipped with a flow-type bactericidal UV cleaner–recirculator AR, which provides constant decontamination inside the cabinet during operation. They are recommended for operations with DNA/RNA amplicons.

UV cleaner–recirculator AR consists of a UV lamp, a fan and dust filters organized in a special body so that a user working with a UV-cleaner box is protected against UV light. Recirculator increases the maximum density of UV light making it sufficiently effective for DNA/RNA inactivation. The UV–recirculator processes 100 UV-cleaner box volumes per hour, creating permanent aseptic conditions of operation inside the UV-cleaner box.

Specially assigned moving tables **T-4** (with wheel locks) with a drawer are available on request.

Advantages of Biosan UV-cleaner boxes:

- Ozone free high density UV decontamination
- Long living UV lamps (9,000 hours average)
- Automatic switch off of UV-lamps when the protective screen is opened
- Bactericidal flow-type recirculator providing permanent decontamination inside UV–cabinet during operation
- Shockproof glass walls
- Low noise, low energy consumption
- Tables for installation of UV-cleaner boxes
- UV-cleaner boxes with the bactericidal UV cleaner–recirculator AR is the patented Biosan solution




CAT. NUMBER

With built in socket and inlet for power cords	With built in socket and inlet for power cords
BS-040109-A06	100-240VAC 50/60Hz Euro socket, Euro plug
BS-040109-A07	100-240VAC 50/60Hz UK socket, UK plug
BS-040109-A08	100-240VAC 50/60Hz US socket, US plug
BS-040109-A24	100-240VAC 50/60Hz AU socket, AU plug
BS-040109-AK	IQ OQ document
BS-040109-BK	PQ document


SPECIFICATIONS

Walls material	Rear: stainless steel, Sides: steel with chemical resistant powder coating Front: glass (EUROGLASS, Germany)
Working surface material	Stainless steel
Open UV-lamp	1 x 25W built-in bactericidal UV-C, TUV 25W 1SL/25
UV radiation level	18 mW / cm ² / sec
Radiation type	UV ($\lambda = 253.7$ nm), ozone-free
Digital time setting of direct UV exposure	1 min–24 hrs / non–stop (increment 1 min)
UV-recirculator	1x25W (efficiency >99% per 1 hour)
Daylight lamp (for working area illumination)	1 x TLD-15W
Thickness of side panels	2 mm
Thickness of upper front panel	8 mm
Thickness of screen	4 mm
Optical transmission	95%
UV protection	>96%, UV-protection film
Working area	645 × 490 mm
Opening dimensions (fully raised protective screen)	645 × 190 mm
Safety features	Automatic open UV-lamp switch off when screen is open
Power outlets inside the unit	Inlet for power cords and 1x Built-in socket, max. 1000 W
Overall dimensions (W×D×H)	700 × 580 × 555 mm
Weight (net / gross)	31.2 / 42 kg
Power consumption	67 W
Nominal operating voltage	100-240 V, 50/60Hz
Optional table	T-4 (800 × 600 × 750 mm)

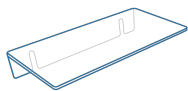

F-1
BS-040104-DK
shelf

F-1, shelf

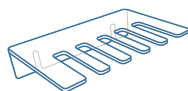

P-5
BS-040104-CK
shelf for pipettes

P-5, shelf for pipettes



T-4
BS-040101-BK
Table

New modular design of
laboratory furniture provides
flexibility and ease of use.



LF-1
BS-050101-BK
Drawer unit

New modular design of
laboratory furniture provides
flexibility and ease of use.



PDS-250
BS-040107-DK
DNA/RNA Decontamination
Solution

DNA/RNA Decontamination
Solution.



PDS-10L
BS-040107-FK
DNA/RNA Decontamination
Solution 10l

DNA/RNA Decontamination
Solution.